




## Dexcom G6 Phone Set-up

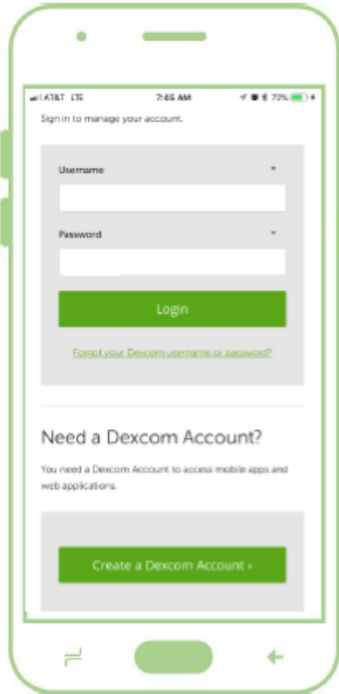

**Step 1:** Log-in to G6 app by entering user name and password provided on the Remote Monitoring Site Worksheet

**Step 2:** Once you have logged in, the app will walk you through a setup wizard

- This includes:
  - Informational videos and screens
  - Safety screens
  - Setting low and high alert levels
  - Entering sensor code
  - Smart device settings (e.g. make sure Bluetooth is on)
  - Entering transmitter serial number (SN)
  - Inserting sensor/attaching transmitter


### Program the following alerts:

- Low alert: 80 mg/dL
- High alert: 300 mg/dL


---

### 4-digit Sensor Code:

- Located on the sensor's adhesive paper
- Unique to that sensor

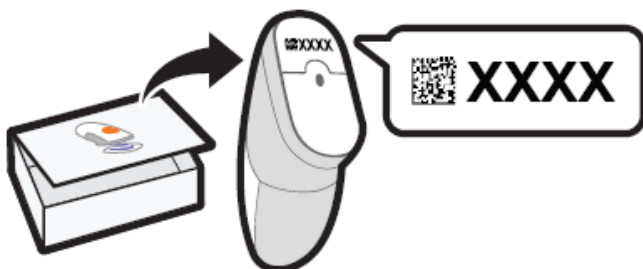

### Transmitter Serial Number (SN)

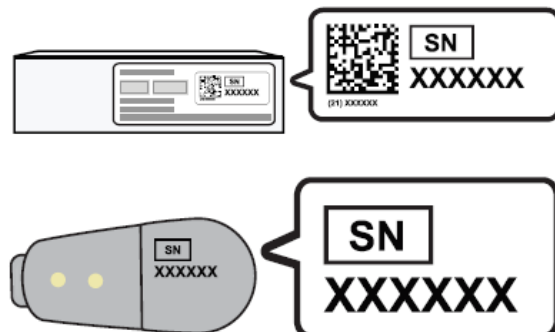

Transmitter SN is on the box and on the bottom of the transmitter—look for **SN**